

Abbreviations and notes:

Col: Colony

Ver: Verified

1:1000 atc corresponds to 250 ng/ μ l

2log is used as a ladder for all gel electrophoresis experiments.

24.06.17

- For pza-HxlR-sfGFP and psb1c3-HxlR PCR, Gibson ve transformation are performed.
- Cultures are prepared.

25.06.17

- psb1c3-hxlR plate contains no colony.
- pza-hxlR-sfgfp plate (plate 1) has 5 colonies.

28.06.17

- Cultures are prepared from colony 1, 2, 3 and 4 of plate 1 (pza-hxlR-sfgfp).


29.06.07

- To the colony 1, 2, 3 and 4 miniprep is performed (pza-hxlR-sfgfp, plate 1).

Col 1- 236.5 ng/ μ l col 2- 57.0 ng/ μ l col 3- 226.6 ng/ μ l col 4- 166.5 ng/ μ l

- The plasmids isolated from miniprep experiment, digestions are prepared with Sall ve BgIII.

Figure 1 : Pza-hxlR-sfgfp, colony(1-4), plate 1, digest Sall + BgIII :

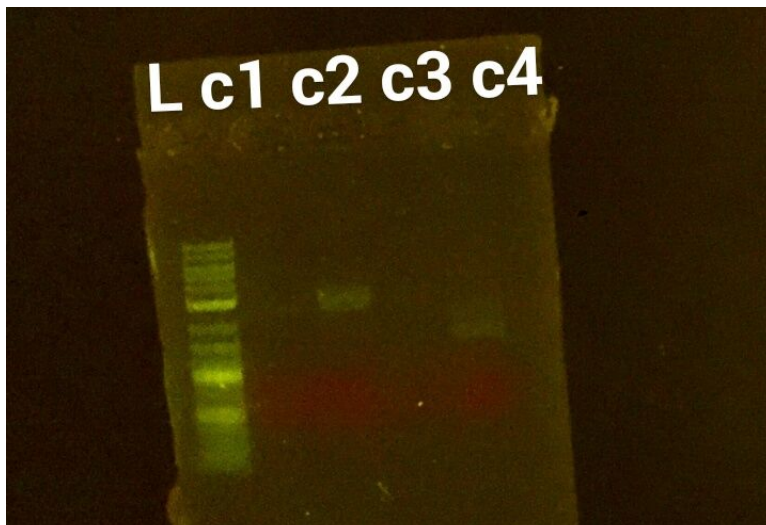

30.06.17

-Miniprep samples from colony 1 and 2 are sent for sequencing (pza-hxIR-sfgfp plate 1).

03.07.17

-Cultures are prepared again from colony 1 and 2 from plate 1.

04.07.17

-Induction is performed for colony 1 and 2. (unind / atc / atc+voc / voc)

1:50 dilution for fresh LB

1:1000 dilution for atc (2-2.5 hours of incubation)

05.07.17

-GFP signals are measured. Unind / atc / voc / atc+voc / pro (c1-5 col1, d1-5 col2)

06.07.17

-Induction is performed for colony 1.

1:1000 atc - 0 / 0.05 / 0.1 / 0.2 / 1 mM formaldehyde

10.07.17

Figure 2: Induction results of pza-hxlR-colony 2. M5 order: LB unind atc 0.05 0.1 0.2 1 (in mM), Induction results of pza-hxlR-colony 2. X axis indicates the GFP measurement after voc addition: 2, 3, 4 and 24 hours. Y axis indicates the normalized GFP reading.

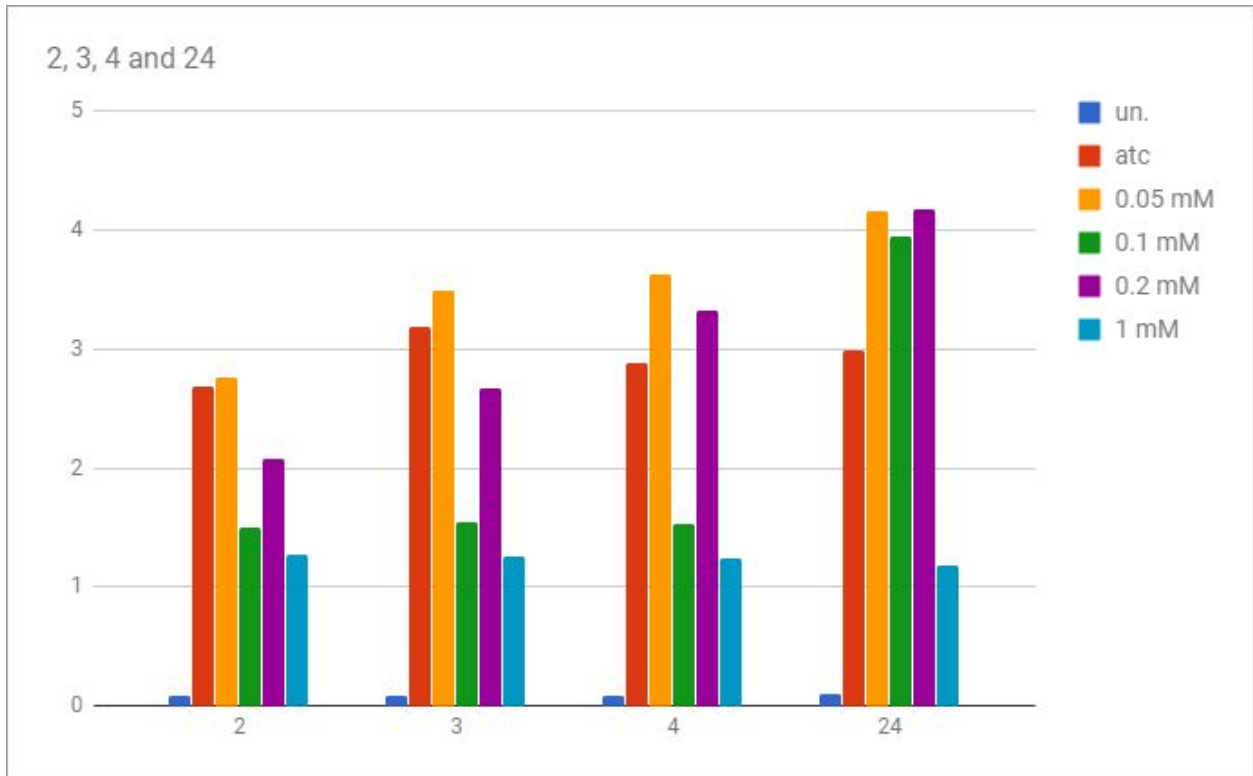


Table 1: Normalized GFP signal data from M5 reader for colony 2 induction

	un.	atc	0.05 mM	0.1 mM	0.2 mM	1 mM
2 h	0.08	2.686	2.764	1.4983	2.0749	1.2638
3 h	0.0815	3.181	3.4851	1.5375	2.6687	1.2503
4 h	0.0886	2.8799	3.6198	1.5309	3.326	1.2403
24 h	0.1048	2.9891	4.1495	3.9383	4.1727	1.184

-Pza-hxlR-sfgfp colony 2 (plate 1) is prepared for cell culture for induction.

11.07.17

-Pza-hxIR-sfgfp-col2 2nd induction: unind. - Atc- 0.05 mM - 0.1 mM - 0.4 mM - 1 mM and pro

Dilution:

For colony 2, 30 ml fresh LB + 600 μ L o.n col2 + 30 μ L cmr

For pro 8 ml fresh LB + 160 μ L + 8 μ L spectinomycin

After 2 hours and 15 minutes incubation 0.38 OD level is achieved.

After 3 hours 25 minutes atc is added.

3 hours later formaldehyde is added. Each sample contains 4 mL of culture.

13 M formaldehyde is diluted to 20 mM. (1ml formaldehyde + 650 μ L ddH₂O)

For 0.05 mM 10 μ L diluted formaldehyde is used.

For 0.1 mM 20 μ L diluted formaldehyde is used.

For 0.4 mM 80 μ L diluted formaldehyde is used.

For 1 mM 200 μ L diluted formaldehyde is used.

12.07.17

Figure 3: Induction results of Col2-pza-hxIR-sfgfp plate 1. GFP measurement order: LB, pro, unind, atc, 0.05, 0.1, 0.4, 1 (mM). X axis indicates the GFP measurement after voc addition: 1, 2, 3, 4 and 24 hours. Y axis indicates the normalized GFP reading.

1st, 2nd, 3rd, 4th and 24th


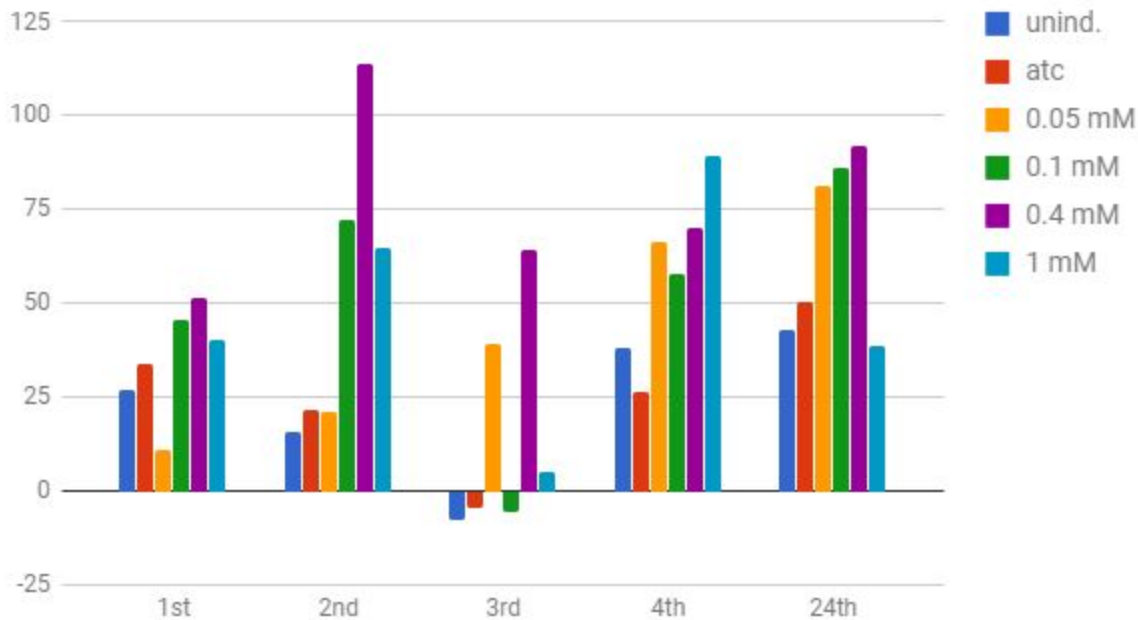


Table 2: Normalized GFP signal data from M5 reader for colony 2 induction

	unind.	atc	0.05 mM	0.1 mM	0.4 mM	1 mM
1st	26.9191	33.601	10.7735	45.3964	51.3835	40.2459
2nd	15.6303	21.6393	21.2892	72.3037	113.9423	64.7239
3rd	-7.813	-4.7165	39.3812	-5.7551	64.3198	4.9212
4th	37.9545	26.1423	66.3795	57.971	69.8467	89.229
24th	42.8865	50.1338	81.3497	86.1527	91.7562	38.5177

Figure 4: XbaI digest for colony 1 and 2 (miniprep products on 29.06.17). On 29.06.17, miniprep products of colony 2 has been verified.

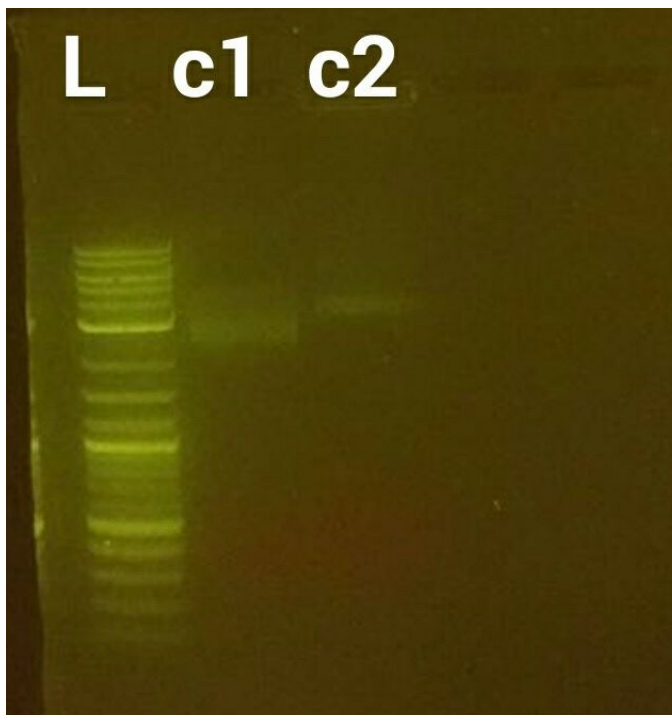

13.07.17

-Miniprep is performed to colony 2. Concentration is 25.7 ng/ μ L

14.07.17

- Transformation is performed to miniprep products from 29.06.17 of colony 2. Since there was little left, it has been vortexed with 10 μ L ddH₂O and used after this. Plates were put to the incubator at 17:00.

Plate name: hxIR verified

15.07.17

Figure 5: Control plate for 14.07.17 transformation


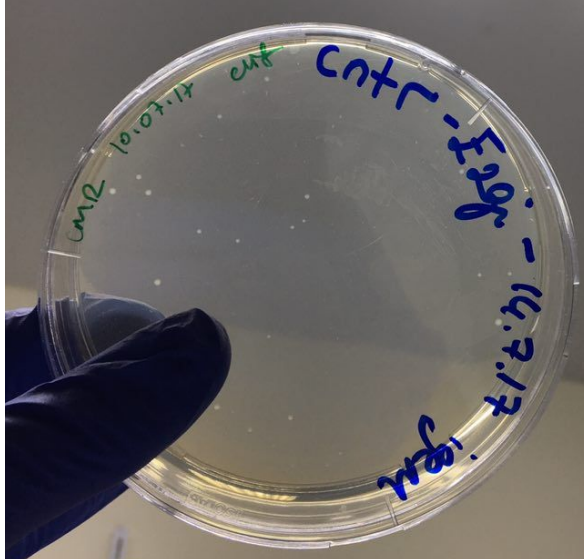


Figure 6: Pza-hxIR-sfgfp (verified)

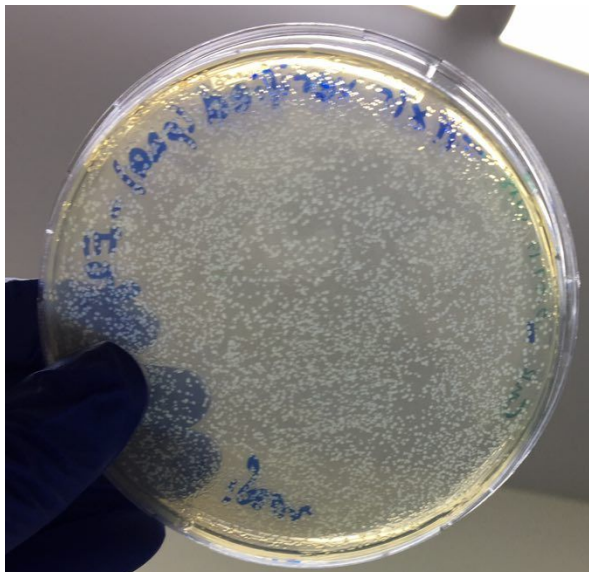

16.07.17

-HxIR ver. col1, ver. col2, ver. col3 and pro are prepared for cell culture. (at 16.30) each 3 mL

17.07.17

-HxIR verified col1-2-3, 2 for each, glycerol has been taken from stock. (200 µL o.n + 200 µL %50 glycerol) One of them is sterile and covered with parafilm. Other is for cell culturing.

-Induction is performed to HxIR ver. Col1, ver. col2 and ver. col3.

For verified col1-2-3 fresh LB:

30 ml LB + 30 µL cmr + 600 µL o.n

Pro:

8 ml LB + 8 µL spek + 160 µL o.n

-For verified col1, col2, col3 and pro amounts of atc:

4 ml - Unind. ---- 0

4 ml - 1:0 ----- 0

4 ml - 1:1,000 --- 4 µL

4 ml - 1:2,500 --- 1.6 µL

4 ml - 1:5,000 --- 0.8 µL

4 ml - 1:10,000 -- 0.4 µL

4 ml - 1:25,000 -- 0.16 µL

4 ml - pro ----- 0

- Formaldehyde amounts for verified col1, col2, col3 and pro: (0.4 mM each induced one)

Dilution: (13 M)(2 µL) = (0.02 M)(1300 µL)

(20 mM)(40 µL) = (0.4 mM)(2000 µL)

4 ml - Unind. ---- 0

2 ml - 1:0 ----- 40 µL

2 ml - 1:1,000 --- 40 µL

2 ml - 1:2,500 --- 40 µL

2 ml - 1:5,000 --- 40 µL

2 ml - 1:10,000 -- 40 µL

2 ml - 1:25,000 -- 40 µL

2 ml - 1:0 ----- 0

2 ml - 1:1,000 --- 0

2 ml - 1:2,500 --- 0

2 ml - 1:5,000 --- 0

2 ml - 1:10,000 -- 0


2 ml - 1:25,000 -- 0

4 ml - pro ----- 0

In order to put HxIR into psB1C3; AlkR, hxIR and psB1C3 PCR was conducted.

AlkR-> Tm: 64 Hxlr-> TM: 64 psb1c3-> Tm:71 Q5


Figure 7: psB1c3 PCR results. Left:AlkR right:Hxlr bottom: pSB1C3 expected band lengths respectively: ~1300 bp,700bp, 2000bp

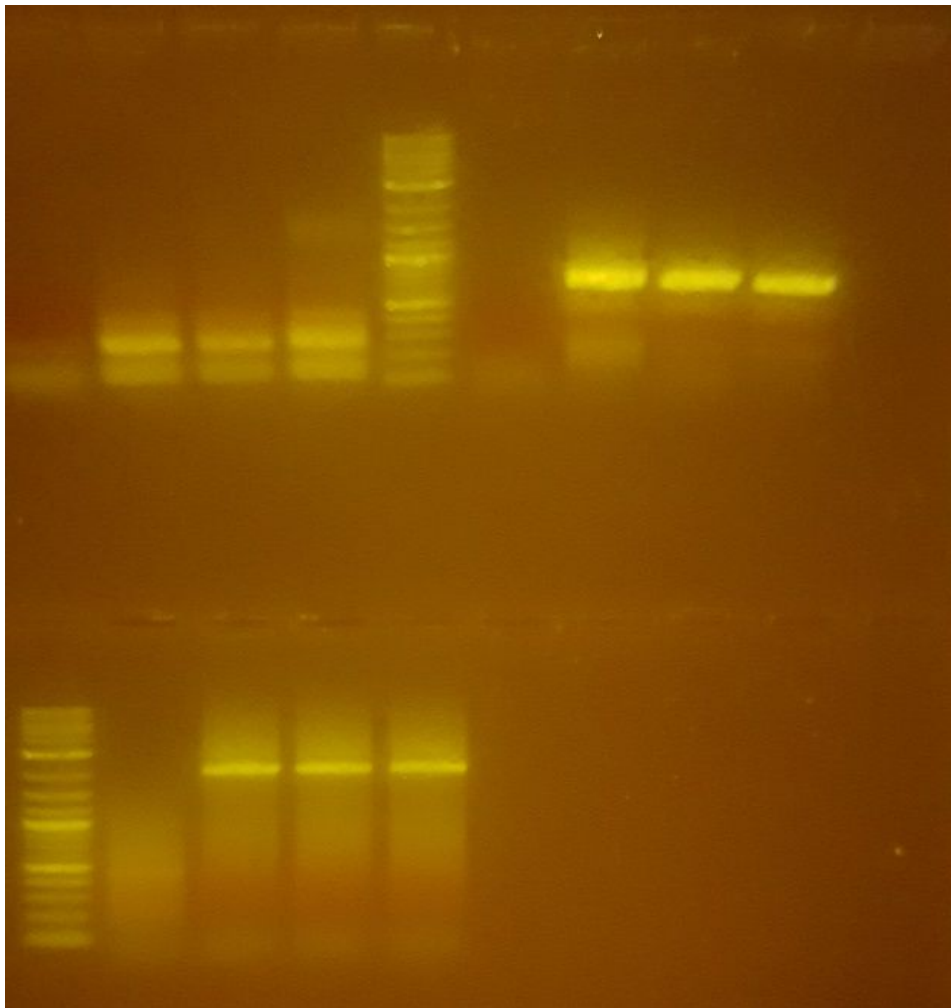

HxIR and psB1C3 were combined by Gibson assay. (3:1, 50 ng:49.42 ng)
19.07.17

For HxIR pSB1C3, 4 colony was selected and colony PCR was performed

Figure 8: Colony PCR results for HxIR psb1c3. Gel: L C 1 2 3 4 3rd and 4th colonies were seen to be have the expected band lengths. (682bp)

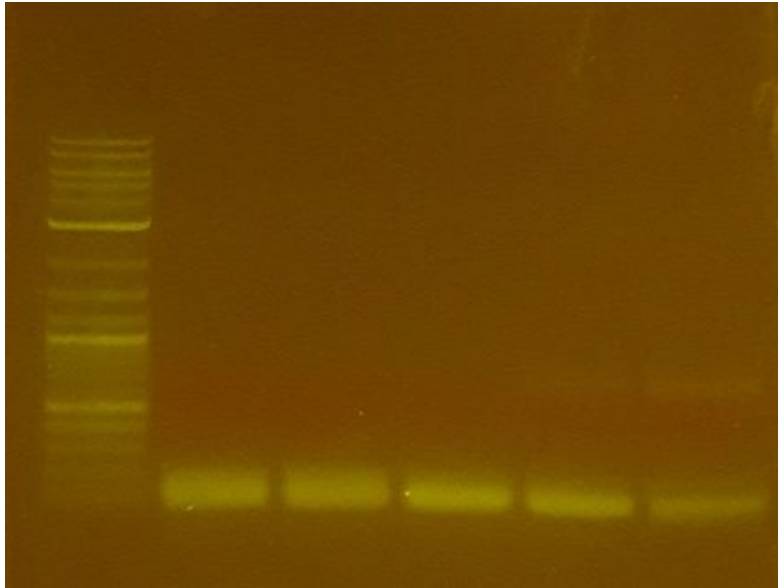


4 HxIR pSB1C3 colony cell culture were prepared.

20.7.17

4 HxIR pSB1C3 miniprep was conducted.

Figure 9: 4 HxIR pSB1C3 and mimR pSB1C3 and AlkR's PCR product was cut with XbaI and SpeI enzymes. Gel: AlkR PCR, HxIR pSB1C3 1-2-3-4, mimR pSB1C3 was all cut by XbaI and SpeI.

1282 bp/2062-496-144 bp/2062 bp-1790 bp

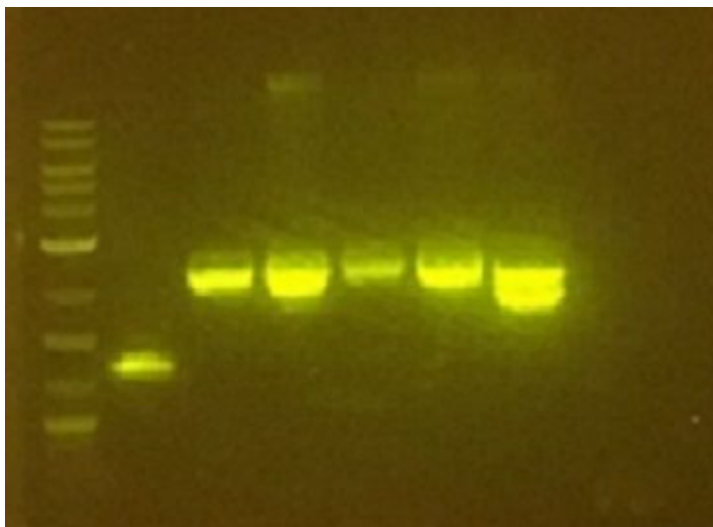


-HxIR pSB1C3 was not verified by the cut with XbaI SpeI enzymes. Cell culture was prepared with newly selected 3 colonies from MimR pSB1C3.

21.07.17

HxIR pSB1C3 5-6-7 miniprep was conducted.

Figure 10: The HxIR's PCR and mimR pSB1C3 and Hxlr pSB1C3 5-6-7 which we already have was digested with XbaI and PstI enzymes. Gel: HxIR psB1C3 5-6-7 , MimR pSB1C3 and HxIR PCR was all digested with XbaI and PstI enzymes.

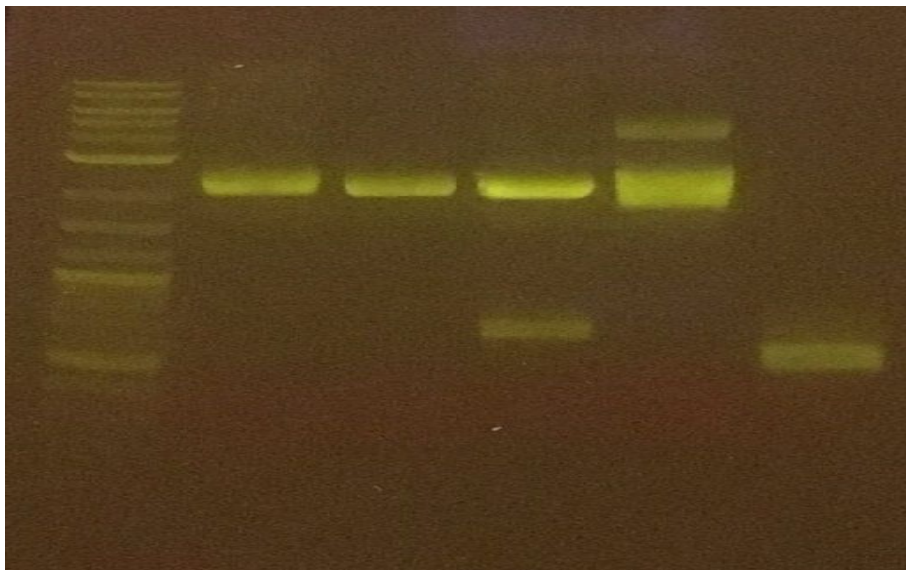

2044-496-162 bp/ 2044-1808 bp/ 658 bp expected bands were not seen. The very end band was seen correct but when we looked zoomed 2 bands were seen. pSB1C3 mimr 5th colony was verified with the sequencing.

24.7.2017

Figure 11: Alkr psb1c3 PCR Control-1-2-3 Hxlr psb1c3 PCR Control- 1-2-3 to be cloned into the pSB1C3.

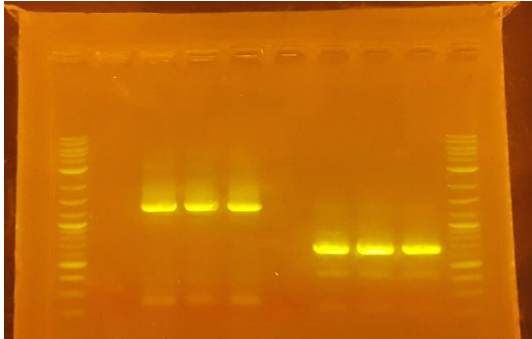


Figure 12: HxIR PCR and AlKR's PCR was digested with XbaI and pstI. But we learned that xbaI also cuts HxIR from the very middle.

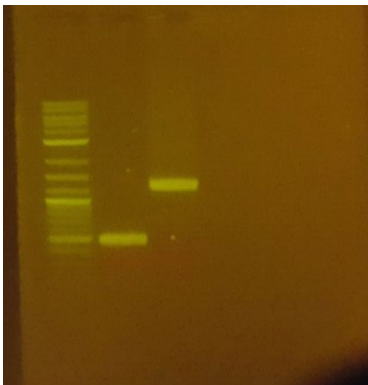


Figure 13: Ver. col1 GFP measurement results. X axis indicates the GFP measurement after voc addition: 1, 2, 3, 4 and 24 hours. Y axis indicates the normalized GFP reading.

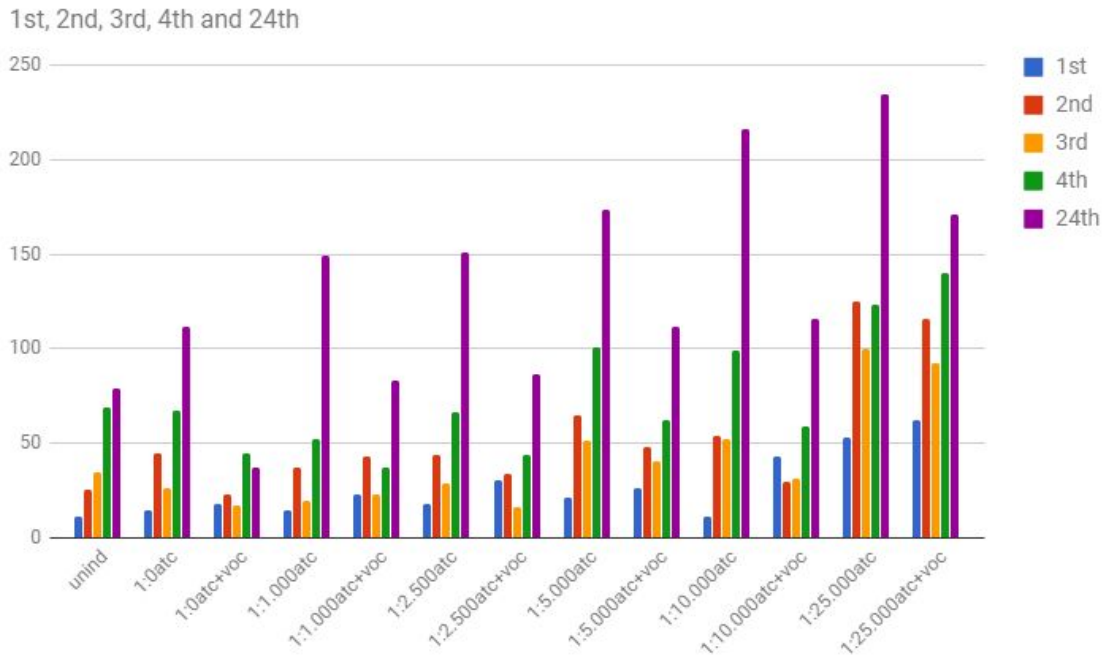


Figure 14: Ver. col2 GFP measurement results. X axis indicates the GFP measurement after voc addition: 1, 2, 3, 4 and 24 hours. Y axis indicates the normalized GFP reading.


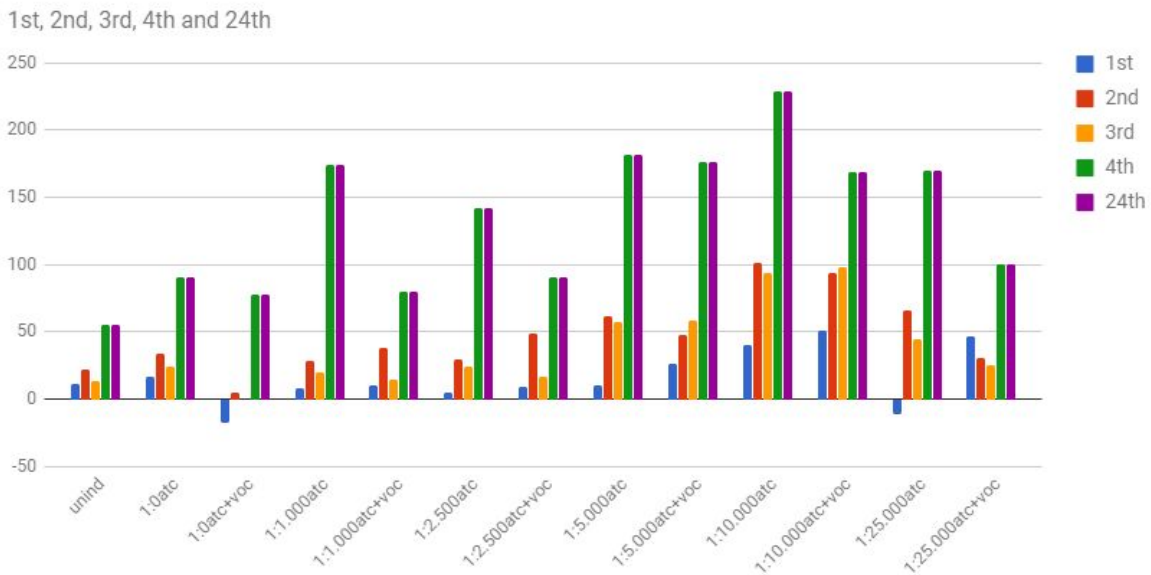


Figure 15: Ver. col3 GFP measurement results. X axis indicates the GFP measurement after voc addition: 1, 2, 3, 4 and 24 hours. Y axis indicates the normalized GFP reading.

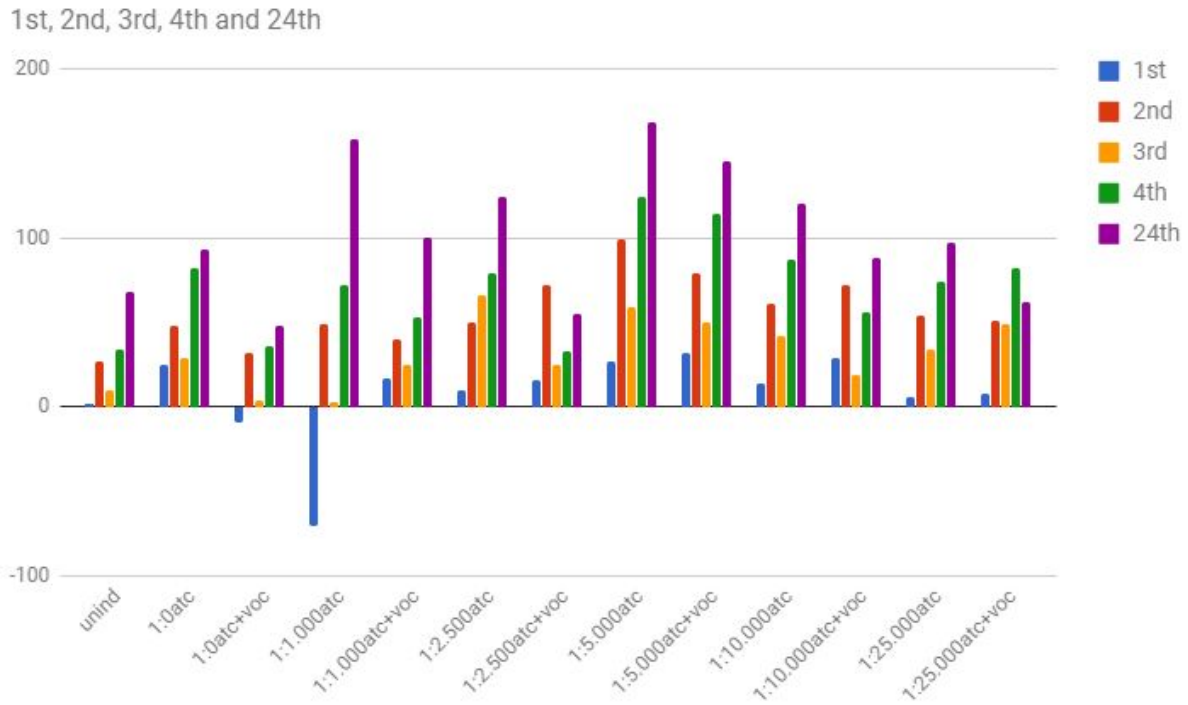

26.7.17

Miniprep was conducted to pSB1C3 HxIR colonies 7-8-9 and digested with EcoRI and PstI. Then its ligation and transformation is completed.

18.8.17

miniprep was conducted to pza HxIR colonies 1-2-3-4. These samples were sent to sequencing.

21.8.17

Histag HxIR digest was conducted with XhoI.

12.10.17


-New induction trials for HxlR verified colony 2:

Table 3 : Experimental design for induction, with 3 replicas and formaldehyde (voc) concentration 1 mM

Hxlr X3	No atc	No atc + voc	1:1000 atc	1:1000 atc + voc	1:10,000 atc	1:10,000 atc + voc	1:25,000 atc	1:25,000 atc + voc	1:50,000 atc	1:50,000 atc + voc
Pro X3	No atc	No atc + voc	1:1000 atc	1:1000 atc + voc	1:10,000 atc	1:10,000 atc + voc	1:25,000 atc	1:25,000 atc + voc	1:50,000 atc	1:50,000 atc + voc

Figure 16: GFP measurement results. X axis indicates samples with different formaldehyde concentrations shown in table 3 and incubation periods: 4 and 16 hours. Y axis indicates the normalized GFP reading.

Note: GFP signal is measured in PBS rather than LB.

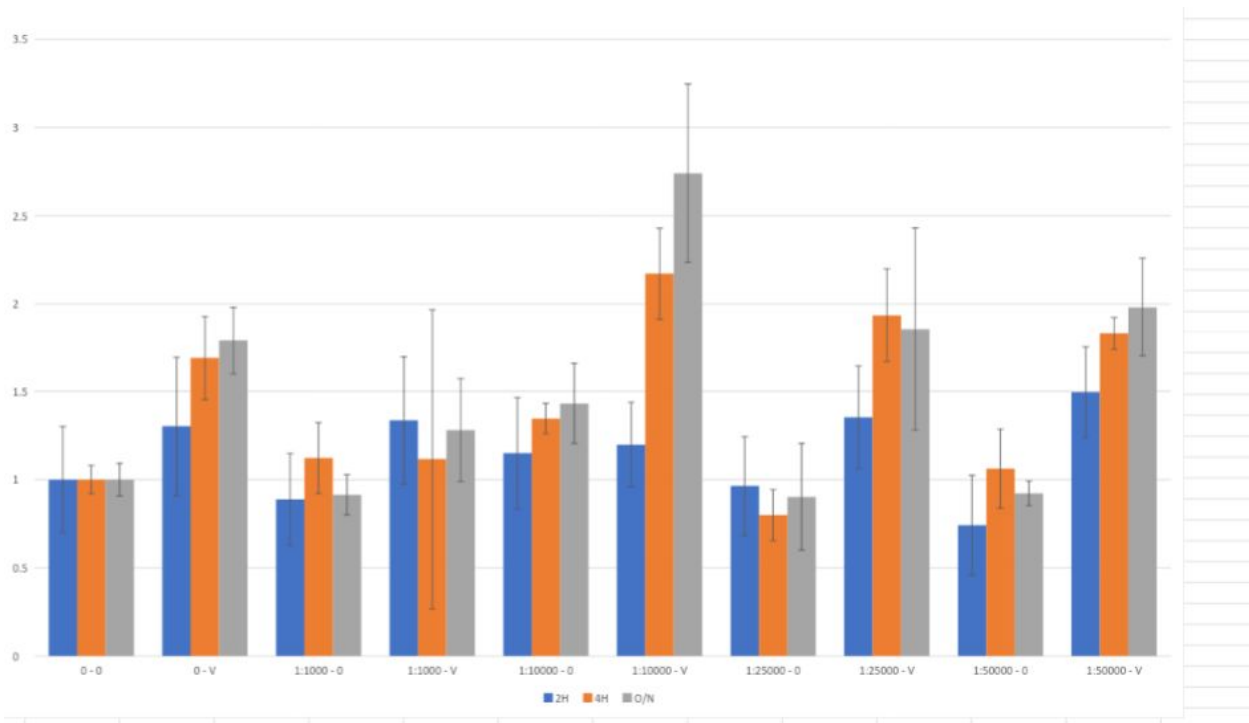

14.10.17

-New induction trials for HxIR verified colony 2:

Table 4 : Experimental desing for induction of HxIR verified col 2 . Gradient formaldehyde (voc) with 3 replicas, 1:10,000 atc

Hxlr X3	atc+ 10 mM voc	atc+ 3.2 mM voc	atc+ 1 mM voc	atc+ 0.32 mM voc	atc+ 0.1 mM voc	atc+ 0.01 mM voc	atc+ 0.001 mM voc	atc+ 0 mM voc
Pro X3	atc+ 10 mM voc	atc+ 3.2 mM voc	atc+ 1 mM voc	atc+ 0.32 mM voc	atc+ 0.1 mM voc	atc+ 0.01 mM voc	atc+ 0.001 mM voc	atc+ 0 mM voc

Figure 17: GFP measurement results after 16 incubation...X axis indicates samples with different formaldehyde concentrations shown in table 4. Y axis indicates the normalized GFP reading.

Note: GFP signal is measured in PBS rather than LB.

